


The ConnectingOntario ClinicalViewer Diagnostic Imaging (DI) portlet provides access to diagnostic images and reports from the Diagnostic Imaging Common Service (DI CS). The condensed portlet, the full-sized portlet and the near full-size portlet display the same columns of information.

Viewing Information about Diagnostic Imaging Reports and Images

- In the Patient Care Views navigation bar:
 - Select **Summary View**  to view the condensed DI portlet
 - Select **Diagnostic Imaging**  or **Summary List View**  to display the full-sized portlet, OR
 - In the **DI portlet**, select **Portlet Expander**  to display the near full-size portlet
- As needed, change the timeline to view results by selecting **Today**, **7D**, **30D**, **3M**, **6M**, **1Y** or **Custom**. (This changes the timeline for all portlets).

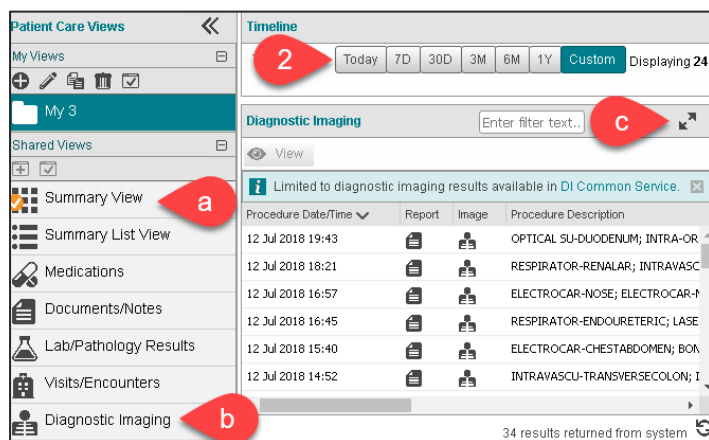

To view additional information about a study:

- By default, the portlet sorts in descending order by **procedure date/time**. Select the **down arrow** to change the order to Ascending. Select the **arrow** again to return to Descending order
- One or more procedures may be associated with an accession. Each row pertains to one accession
- Select **Report**  to view the report in the Document Viewer
- Select **Image**  to view the related diagnostic imaging study/studies in the Image Viewer
- Select **Print**  to print the contents of the DI portlet

Diagnostic Imaging							
View							
Procedure Date/Time	Image	Procedure Description	Additional Information	Report Status	Organization	Accession Number	
12 Jul 2018 08:10		POSITRON E-PROFFEMVL		Final	North York Gene	ACC IM0072	
12 Jul 2018 07:15		DIGITAL RA-PELVIS		Final	North York General Hospital	ACC IM0069	
12 Jul 2018 06:05		OPTICAL SU-BUTTOCKR; OPTICAL SU-BUTTOCKR; O...		Amended (A)	North York General Hospital	ACC IM0074	
12 Jul 2018 03:48		RA-SIGMOID; BONE DENS-ENDORECTAL; B...	(Hypopharynx, Bone Densitometry (ultras...	Final	North York General Hospital	ACC IM0054	
12 Jul 2018 03:45		ULTRASOUND-ENDOVASCULAR; INTRAVASCU-CTSPI...	(Intravascular Optical Coherence Tomogr...	Final	North York General Hospital	ACC IM0079	

- Hover over **Procedure Description** to view the procedure names for each report or study associated with the accession. Each bullet in the tooltip is associated with one report or study
- Additional information** such as modality, anatomic region and laterality display when provided by the source system

Diagnostic Imaging							
View							
Procedure Date/Time	Report	Image	Procedure Description	Additional Information	Report Status	Organization	Accession Number
12 Jul 2018 01:28			ELECTROCAR-NOSE; THERMOGRAP-SMALLINTESTIN...	(Thermography, Perineum); (Thermograp...	Final	North York General Hospital	ACC IM0057
12 Jul 2018 00:58			ELECTROCAR-NOSE • THERMOGRAP-SMALLINTESTINE • THERMOGRAP-SMALLINTESTINE • THERMOGRAP-SMALLINTESTINE	(Magnet • (Thermography, Perineum) • (Thermography, Perineum) • (Perineum, Thermography)		North York General Hospital	ACC IM0080


Viewing Reports in the Document Viewer

To view reports:

1. Select **Report**  in the appropriate row
The Document Viewer opens in a separate window

Report Tips:

2. Each report is visible in its own tab
3. Amended DI reports display with a **red** reporting status label and the Amended label. Amended reports replace the original report, although both display in the Document Viewer
4. Print the report using **Print Current Tab**
5. Close the report window using **X** or **Close**

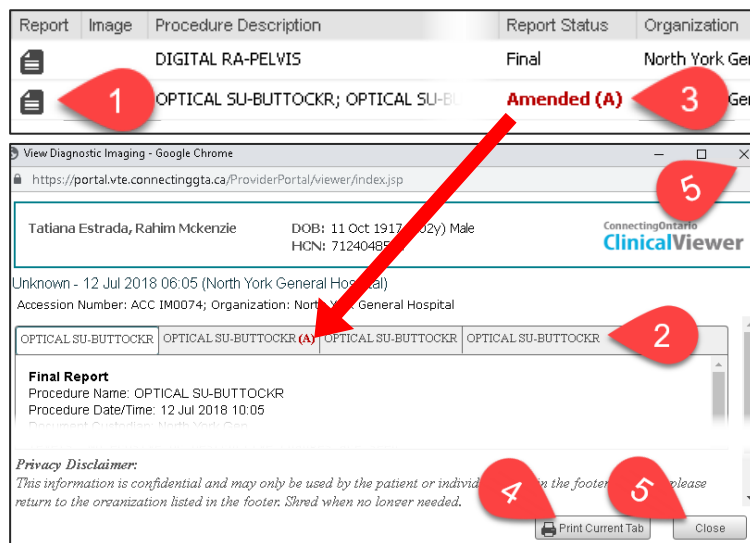

Viewing Studies in the Image Viewer

Note: Your desktop or laptop monitor may not produce the same viewing quality as a PACS (Picture Archiving and Communication System) Viewer

To view diagnostic imaging studies:

1. Select **Image**  in the appropriate row
The Image Viewer opens in a separate browser tab and automatically logs you in

Image Viewer Tips:



2. The **red border** identifies the active image/study
3. **Navigate** to the next series by page within the study group
4. Display thumbnails of all available series in the accession, allowing you to scroll through all

Note: For detailed information on the use of the Image Viewer, select **Help** on the toolbar

Displaying Images from Multiple Accessions Simultaneously

Images for up to four accessions can be viewed at the same time:

1. Select **Image** on the row of the first accession to be viewed

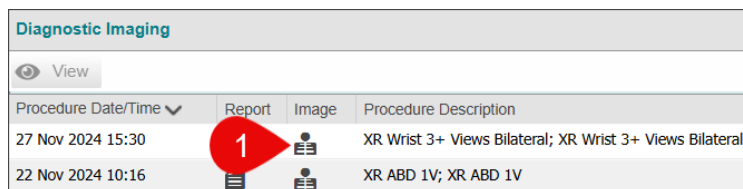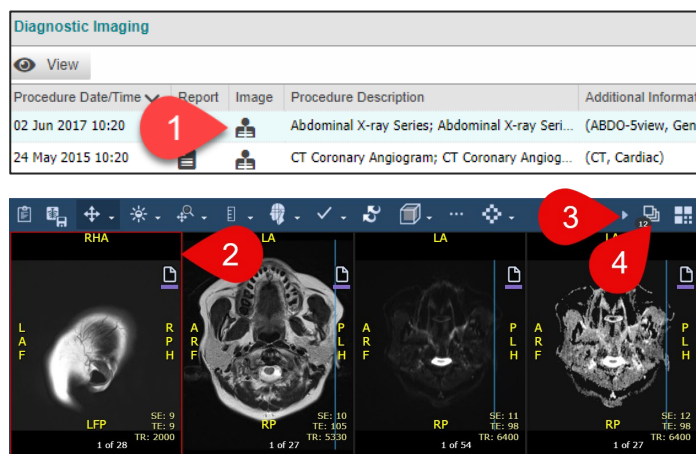


How to Use the Diagnostic Imaging Portlet

2. On the Image Viewer toolbar, select **Layout** and select the desired layout from the list
 - The Image Viewer window displays the number of viewports selected with the current study in the active one, indicated by a red border




3. Select the viewport the next study should display in. On the Image Viewer toolbar, select **Clipboard**
 - The Studies menu displays, showing all accessions containing images for the patient


Note: The studies/accessions listed in the Studies menu **is not restricted** to the date range currently selected in the ClinicalViewer timeline.

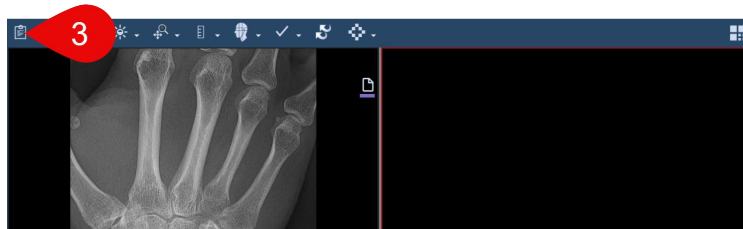

4. Select the accession to display in the active viewport

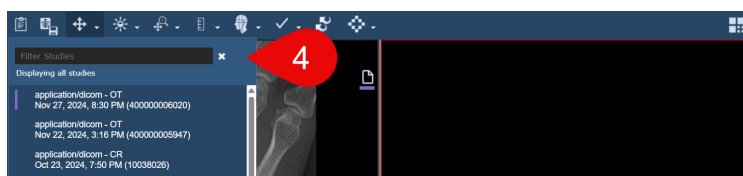

- The study displays in the pane. All tools on the toolbar are available for the active viewport. Activate a different viewport by selecting it.

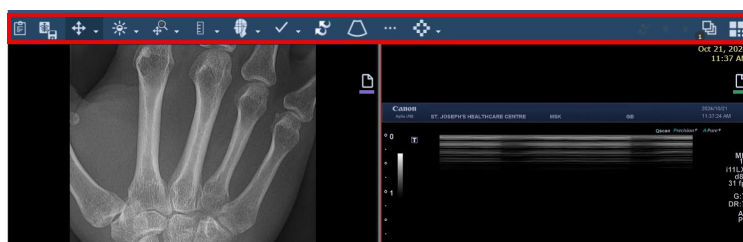


Displaying Images and Reports Simultaneously

To view a report simultaneously on screen with a DI study:


1. Select the **Image** icon in the appropriate row of the DI portlet
2. Return to the ClinicalViewer by selecting its browser tab or window

Procedure Date/Time	Report	Image	Procedure Description
02 Jun 2017 10:20			Abdominal X-ray Series;

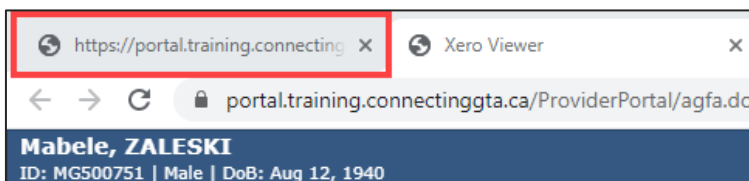


How to Use the Diagnostic Imaging Portlet

3. Select the **Report** icon in the appropriate row within the DI portlet of the ClinicalViewer

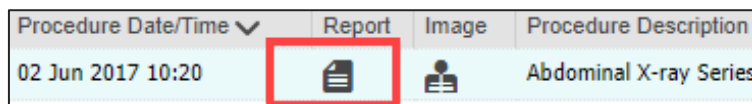


Return to the Image Viewer by selecting its browser tab or window

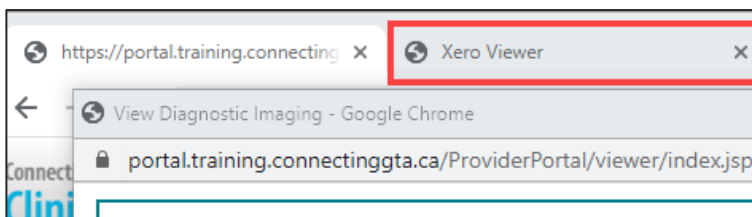


4. Select the report's document viewer by*:
 - a) hovering over the browser icon on the taskbar, and
 - b) selecting the report's window

* The specific steps to follow are dependent on the configuration of your internet browser. If assistance is required, contact your local IT resource.

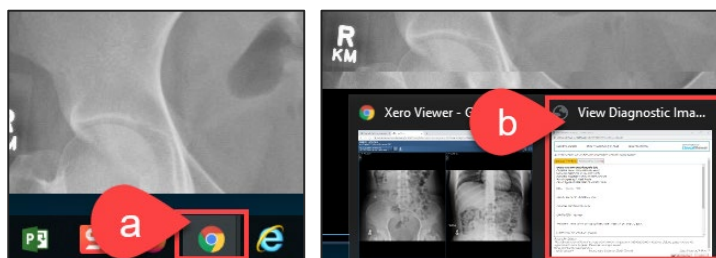

5. Resize and reposition the windows so that both display on screen

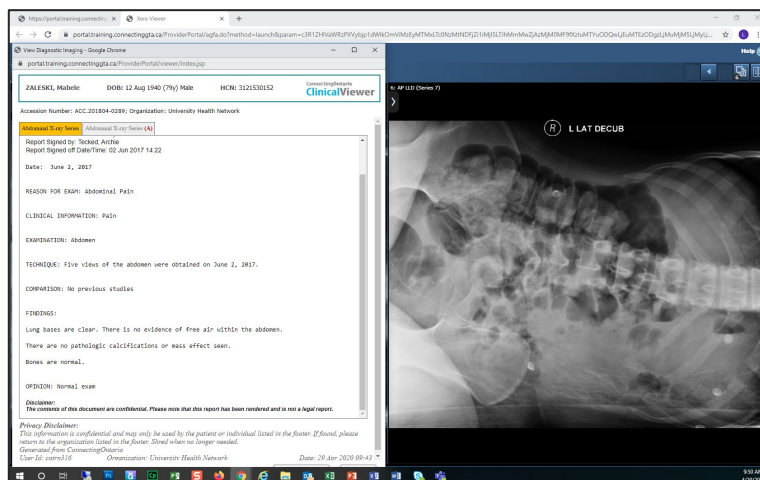

Consent for the Diagnostic Imaging Portlet

- A message displays stating that some or all data in DI CS has been blocked due to the patient's consent directive.
- Currently, consent blocks cannot be overridden in the DI portlet nor can blocked DI information be retrieved from DI CS.



Refer to the [Privacy and Security Information tip sheet](#) for further information about consent directives